





# Characteristic X-rays of tungsten



Physics	Modern Physics	Production & use of X-rays
 Difficulty level hard	 Group size 2	 Preparation time 45+ minutes
		 Execution time 45+ minutes

**PHYWE**  
excellence in science

# General information

## Application

**PHYWE**  
excellence in science

Setup

Most applications of X rays are based on their ability to pass through matter. Since this ability is dependent on the density of the matter, imaging of the interior of objects and even people becomes possible. This has wide usage in fields such as medicine or security.

## Other information (1/2)

**PHYWE**  
excellence in science**Prior****knowledge****Main****principle**

The prior knowledge required for this experiment is found in the Theory section.

An X-ray tube with a tungsten anode generates X-radiation that is selected with the aid of a monocrystal as a function of the Bragg angle. A Geiger-Müller counter tube measures the intensity of the radiation. The glancing angles of the characteristic X-ray lines are then used to determine the energy.

## Other information (2/2)

**PHYWE**  
excellence in science**Learning****objective****Tasks**

The goal of this experiment is to get to investigate the characteristic X-radiation of tungsten.

- Record the X-ray spectrum that is emitted by the tungsten anode as a function of the Bragg angle and with the aid of the LiF monocrystal that is used as the analyser.
- Determine the energy values of the characteristic X-rays of tungsten based on the spectra and compare them with the values that were determined based on the corresponding energy-level diagram.

## Theory (1/3)

When electrons hit the metallic anode of the x-ray tube with a high kinetic energy, x-rays with a continuous energy distribution (bremsstrahlung) are generated. The bremsstrahlung is superimposed by additional discrete lines. If an atom of the anode material is ionised, for example, in a deeper shell by an electron impact, an electron from a higher shell can take the now free place while emitting inter alia an x-ray quantum. The energy of this x-ray quantum corresponds to the energy difference of the two levels that are involved in this process. Since the energy difference is atom-specific, the radiation that is generated by this process is also called the characteristic x-radiation.

When an x-ray with the wavelength  $\lambda$  hits the group of lattice planes of a monocrystal at the glancing angle  $\theta$ , the rays that are reflected by the lattice planes only interfere in a constructive manner, if their path difference  $\Delta$  corresponds to an integer of the wavelength (see Fig. 1). This condition is represented by the so-called Bragg's law:

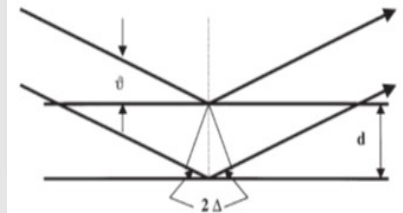


Fig. 1: Bragg scattering on the lattice planes

## Theory (2/3)

$$2d \sin(\theta) = n\lambda \quad (1)$$

(d: interplanar spacing;  $n = 1, 2, 3, \dots$ )

If the interplanar spacing  $d$  is known, the wavelength  $\lambda$  can be determined with the aid of the glancing angle  $\theta$ . The energy of the radiation then results from:

$$E = h \cdot f = \frac{hc}{\lambda} \quad (2)$$

When combining (1) and (2), we obtain:

$$E = \frac{n \cdot h \cdot c}{2d \cdot \sin(\theta)} \quad (3)$$

$$\text{Planck's constant } h = 6.6256 \cdot 10^{-34} \text{ Js}$$

$$\text{Velocity of light } c = 2.9979 \cdot 10^8 \frac{\text{m}}{\text{s}}$$

$$\text{Interplanar spacing LiF (200) } d = 2.014 \cdot 10^{-10} \text{ m}$$

$$\text{Equivalent } 1 \text{ eV} = 1.6021 \cdot 10^{-19} \text{ J}$$

## Theory (3/3)

Figure 2 shows the energy level diagram of the tungsten atom. Since the energy of the K-shell is approximately 70 keV, the maximum available energy of the primary beam of the x-ray unit of 35 keV is insufficient for exciting the K-shell. An ionization is only possible for the L-level. For clarity, Figure 2 shows the L-transitions for the dipole radiation (see 4) that are possible in accordance with the quantummechanical selection rules only up to the N-shell. The quadrupole radiation with a much lower intensity can be neglected.

$$\Delta l = \pm 1 \text{ and } \Delta j = 0, \pm 1 \text{ selection rules for the dipole radiation} \quad (4)$$

( $l$  = orbital angular momentum,  $j$  = total angular momentum)

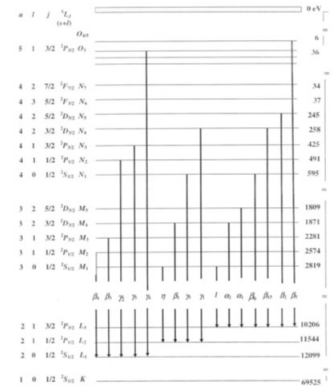


Fig. 2: Schematic representation of the energy level of tungsten (Z = 74)

## Equipment

Position	Material	Item No.	Quantity
1	<a href="#">XR 4.0 expert unit, 35 kV</a>	09057-99	1
2	<a href="#">XR 4.0 X-ray goniometer</a>	09057-10	1
3	<a href="#">XR4 X-ray plug-in W tube</a>	09057-81	1
4	<a href="#">XR 4.0 X-ray solid state physics upgrade set</a>	09125-88	1





# Setup and Procedure

## Setup

Connect the goniometer and the Geiger-Müller counter tube to their respective sockets in the experiment chamber (see the red markings in Fig. 3). The goniometer block with the analyser crystal should be located at the end position on the right-hand side. Fasten the Geiger-Müller counter tube with its holder to the back stop of the guide rails. Do not forget to install the diaphragm in front of the counter tube (see Fig. 4). Insert a diaphragm tube with a diameter of 2 mm into the beam outlet of the tube plug-in unit.

**For calibration:** Make sure, that the correct crystal is entered in the goniometer parameters. Then, select "Menu", "Goniometer", "Autocalibration". The device now determines the optimal positions of the crystal and the goniometer to each other and then the positions of the peaks.

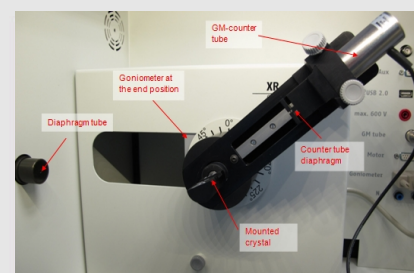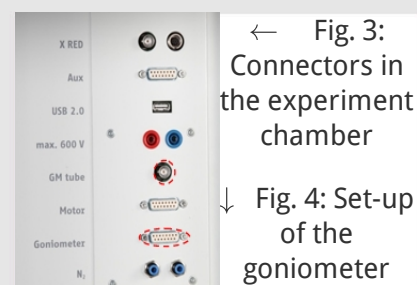

## Procedure (1/3)



- Connect the X-ray unit via the USB cable to the USB port of your computer (the correct port of the X-ray unit is marked in Figure 5).
- Start the “measure” program. A virtual X-ray unit will be displayed on the screen.
- You can control the X-ray unit by clicking the various features on and under the virtual X-ray unit. Alternatively, you can also change the parameters at the real X-ray unit. The program will automatically adopt the settings.



Fig. 5: Connection of the computer

## Procedure (2/3)

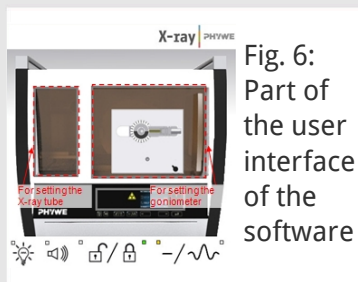


Fig. 6: Part of the user interface of the software

- Click the experiment chamber (see the red marking in Figure 6) to change the parameters for the experiment. Select the parameters as shown in Figure 7 for the LiF crystal.
- If you click the X-ray tube (see the red marking in Figure 6), you can change the voltage and current of the X-ray tube. Select the parameters as shown in Fig. 8.

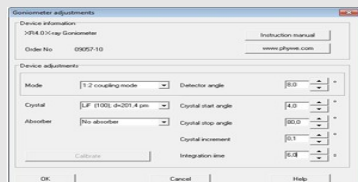


Fig 7: Settings of the goniometer

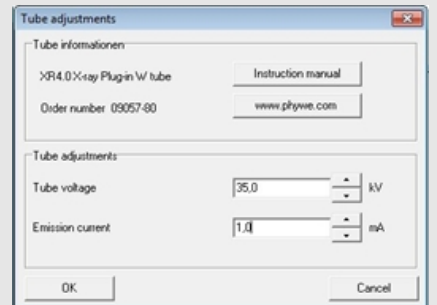



Fig 8: Voltage and current settings



## Procedure (3/3)

**PHYWE**  
excellence in science

- Start the measurement by clicking the red circle:
- After the measurement, the following window appears:

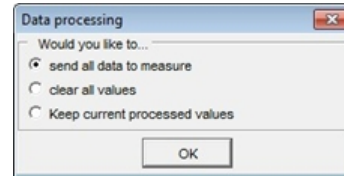

- Select the first item and confirm by clicking OK. The measured values will now be transferred directly to the "measure" software.
- At the end of this manual, you will find a brief introduction to the evaluation of the resulting spectra.

### Overview of the settings of the goniometer and X-ray unit:

- 1:2 coupling mode
- Gate time 5 - 6 s; angle step width 0.1°
- Scanning range 4°-80° (LiF monocrystal)
- Anode voltage  $U_A = 35$  kV; anode current  $I_A = 1$  mA

**PHYWE**  
excellence in science



## Evaluation

## Task 1

**PHYWE**  
excellence in science

### Analyse the intensity of the tungsten X-radiation as a function of the Bragg angle and with the aid of a LiF monocrystal.

Figure 9 shows the X-ray spectrum of tungsten that was analysed with a LiF monocrystal. Well-defined lines are superimposed on the continuous bremsstrahlung spectrum. The glancing angles of these lines are unaffected by the anode voltage. This identifies them as characteristic X-ray lines. The evaluation (see the table) shows that there are only first-order lines in the angular range of  $10^\circ < \theta < 30^\circ$ . They reach the highest intensity level. The range of  $30^\circ < \theta < 80^\circ$  includes the lines with  $n = 2$  and  $n = 3$ . The separation of lines number 2 and 10 into  $\alpha_1$  and  $\alpha_2$  or  $\gamma_2$  and  $\gamma_3$  can only be observed in the range with  $n = 2$ . Line number 11 can be clearly assigned to the  $K_\alpha$  line of copper. The small circular tungsten anode is actually embedded in a cylindrical copper rod that is also partly hit by the electrons.

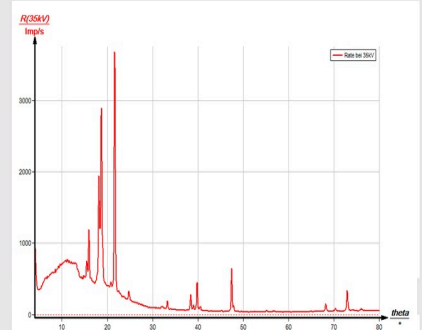


Fig. 9: X-ray spectrum of tungsten

## Task 2

**PHYWE**  
excellence in science

### Determine the energy values of the characteristic X-rays of tungsten based on the spectra and compare them with the values that were determined based on the corresponding energy-level diagram.

Column B of the table shows the glancing angles  $\theta$  that were determined with the aid of Figure 9 and also the energy values for the characteristic X-ray lines of tungsten that were calculated with the aid of equation (4). The wavelength  $\lambda$  and the corresponding energies  $E_{\text{exp}}$  that were calculated with (1) and (3) are shown in columns D and E. Column H shows the energy values  $E_{\text{lit}}$  that were calculated with the aid of Figure 8. The correspondence between the two energy values provides evidence concerning the assignment of the lines to the various transitions. As expected, there are only those lines that comply with the selection rules. It is not possible to discern all of the possible transition, since for some of them the intensity is too low.

## Task 2 (part 2)

A	B	C	D	E	F	G	H	A	B	C	D	E	F	G	H
Line	$\theta$ [°]	$n\lambda$ [pm]	$E_{\text{exp}}$ [eV]	Line	Transition	$E_{\text{lit}}$ [eV]		Line	$\theta$ [°]	$n\lambda$ [pm]	$E_{\text{exp}}$ [eV]	Line	Transition	$E_{\text{lit}}$ [eV]	
1	14.69	1102.15	12138	$\gamma_4$	$L_1O_3$	12063		10	21.47	1147.43	8409	$\alpha_{1/2}$	$L_3M_5/L_3M_4$		
2	15.23	1105.81	11717	$\gamma_{3/2}$	$L_1N_3/L_1N_2$			11	22.51	1154.21	8040	Cu-K	$\alpha_{1/2}$		
3	15.74	1109.27	11346	$\gamma_1$	$L_2N_4$	11286		12	24.57	1167.49	7402	I	$L_3M_1$	7387	
4	16.28	1112.92	10980	$\gamma_5$	$L_2N_1$	10949		13	31.80	2106.13	11682	$\gamma_3$	$L_1N_3$	11674	
5	17.92	1123.94	10003	$\beta_2$	$L_3N_5$	9961		14	32.01	2106.76	11613	$\gamma_2$	$L_2N_4$	11608	
6	18.21	1125.87	9849	$\beta_3$	$L_1M_3$	9818		15	33.03	2109.79		$\gamma_1$	$L_3N_5$	11286	
7	18.47	1127.61	9716	$\beta_1$	$L_2M_4$	9673		16	38.12	2124.33	9972	$\beta_2$	$L_1M_3$	9961	
8	18.79	1129.74	9556	$\beta_4$	$L_1M_2$	9525		17	38.80	2126.20	9824	$\beta_3$	$L_1M_2$	9818	
9	20.60	1141.72	8748	$\eta$	$L_2M_1$	8725		18	39.52	2128.16	9674	$\beta_1$	$L_2M_4$	9673	

## Task 2 (part 3)

A	B	C	D	E	F	G	H
Line	$\theta$ [°]	$n\lambda$ [pm]	$E_{\text{exp}}$ [eV]	Line	Transition	$E_{\text{lit}}$ [eV]	
19	40.24	2130.10	9529	$\beta_4$	$L_1M_2$	9525	
20	47.12	2147.58	8401	$\alpha_1$	$L_3M_5$	8397	
21	47.58	2148.68	8339	$\alpha_2$	$L_3M_4$	8335	
22	54.88	3109.71	11300	$\gamma_1$	$L_2N_4$	11286	
23	56.47	2167.88	7385	I	$L_3M_1$	7387	
24	67.90	3124.28	9976	$\beta_2$	$L_3N_5$	9961	
25	70.09	3126.12	9831	$\beta_3$	$L_1M_3$	9818	
26	72.66	3128.04	9683	$\beta_1$	$L_2M_4$	9673	
27	75.79	3130.03	9535	$\beta_4$	$L_1M_2$	9525	



Table 1

## Note

**PHYWE**  
excellence in science

### "measure" software

With the "measure" software, the peaks in the spectrum can be determined rather easily:

- Click the button  "Mark" and select the area for the peak determination.
- Click the button  "Peak analysis".
- The window "Peak analysis" appears (see Fig. 10). Then, click "Calculate".
- If not all of the desired peaks (or too many of them) are calculated, readjust the error tolerance accordingly.
- Select "Visualise results" in order to display the peak data directly in the spectrum.

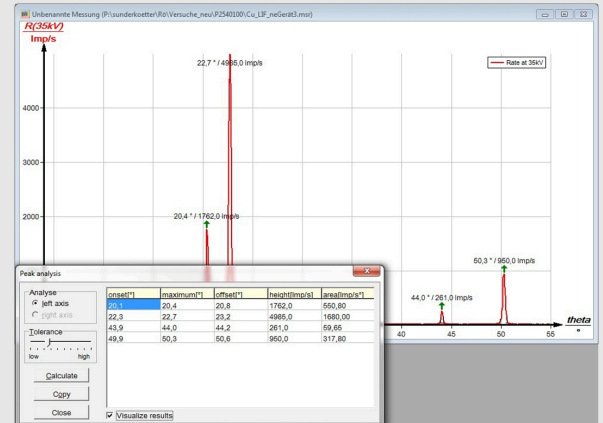


Fig. 10: Automatic peak analysis with "measure"